

4.5 INCH HANDSFREE 2.5X~4X Neck Magnifier with LED Light

Model: ELC-M110LN

INTRODUCTION


This product consists of lighting box, supporter hanging rope and lens. The lighting box contains 2 ultra bright LED lights powered by two included CR2016 button cells. To replace cells, push the cell cover out slightly with finger, install the cells according to the polarity shown on them, and then cover them.

QUICK START

Hang the rope on the neck, put the supporter against the chest, adjust the distance between lens and observed object with the hanging rope in order to obtain clear image. At night or in the dark condition it is necessary to turn on the LED light to get best vision.

CAUTIONS

Don't touch the lens with hands, and always keep the lens clean. If you are not using the magnifier for long time, take the cells out of the magnifier.

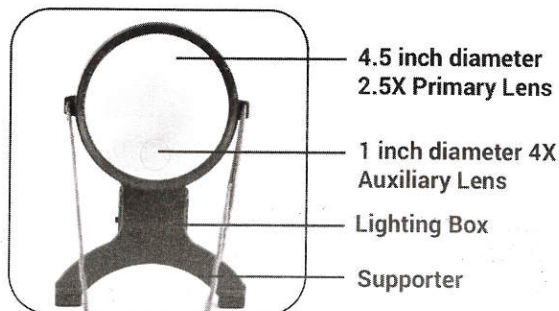

EasyLifeCare®

FEATURES

The chest support LED lighting magnifier can free your two hands so that the observed object can be operated for example during embroidering, knitting, reading and writing. Additionally it is portable and easy to carry. The magnifier lens consists of 2.5X primary lens and 4X auxiliary lens. The product has a function of lighting with high brightness. Two lithium cells (CR2016) can light the two LED bulbs continuously for 20 hours.

CHANGE BATTERIES

1. Hold the body of the magnifier by left hand, clip the position show in the circle in Fig.1 by thumb of the right hand, turn the cell cover and open the cell box following the direction of arrow;
2. Insert two 3V button cells with negative polarity upward in the cell box, as shown in Fig.2;
3. After the button cells are installed, cover the cell box, as Show in Fig.3.

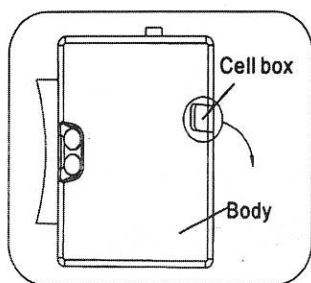


Fig.1

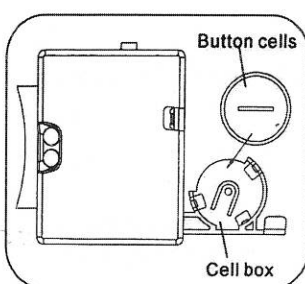


Fig.2

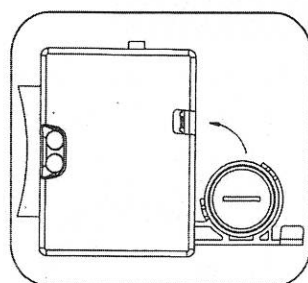


Fig.3

Contact Us : support@easylifecare.com

Made In China